Digital Body Weight Scale

Model: EB9380H

User Manual

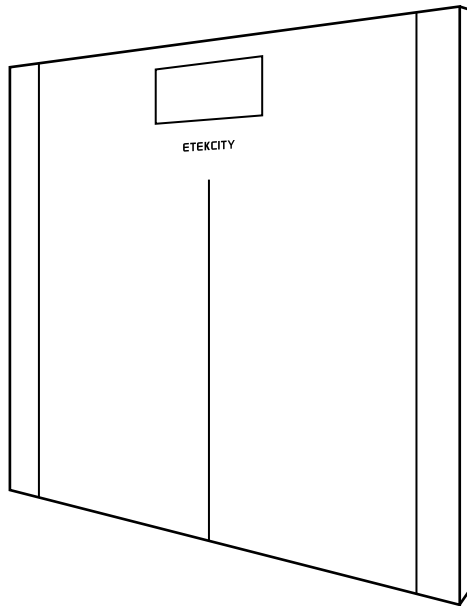
Table of Contents

Package Contents	2
Specifications	2
Safety Information	8
Function Diagram	10
Before First Use	11
Operation	12
• Using the Scale	12
• Body Measuring Tape	13
• Replacing the Batteries	13
• Cleaning the Scale	14
• Calibrating the Scale	14
Troubleshooting	15
Warranty Information	16
Customer Support	17

Package Contents

- 1 x Digital Body Weight Scale
- 1 x Body Measuring Tape
- 2 x 1.5V AAA Batteries (Pre-Installed)
- 1 x Quick Start Guide

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Thank you for purchasing the Digital Body Weight Scale by Etekcity.

If you have any questions or concerns, please reach out to our helpful Customer Support Team at support.eu@etekcity.com. We hope you enjoy your new scale!

BECOME AN ETEKCITIZEN

Get exclusive deals, giveaways, and product registration. Better products for better living.

Find us at etekcity.com

Safety Information

Please read and follow all instructions and safety guidelines in this manual.

General Use

- Check the scale before using. **Do not** use the scale if it is damaged in any way.
- **Do not** strike or bang on the scale and **do not** use the scale if there is a slight crack on the glass. Due to the nature of tempered glass, the glass will shatter if there is a slight crack. If you see a crack, contact [Customer Support](#) (see page 17).
- The glass platform is slippery when wet. **Always** make sure the platform and your feet are dry before stepping onto the scale.
- Place the scale on a hard, flat surface.
- To weigh safely and accurately, stand with your feet shoulder-width apart while on the scale. **Do not** stand on the edge of the scale or jump on the scale.
- Avoid standing on the scale for a long period of time.
- **Do not** overload the scale. Maximum weight capacity is 180 kg / 28 st / 400 lb.
- Keep the scale in a cool area with low humidity. Keep the scale away from heat sources, such as ovens or heaters.
- Store the scale horizontally, not vertically.
- For household use **only**. This scale is not intended for commercial use.
- Handle with care.



Batteries

- **Do not** mix new and old batteries.
- If the scale is not in use for a long period of time, remove the batteries to extend the lifetime of the scale.

SAVE THESE INSTRUCTIONS

Function Diagram

- | | | | |
|----|-------------------|----|-----------------------|
| A. | Platform | D. | Battery Compartment |
| B. | Display | E. | Unit Switching Button |
| C. | Anti-Skid Padding | | |

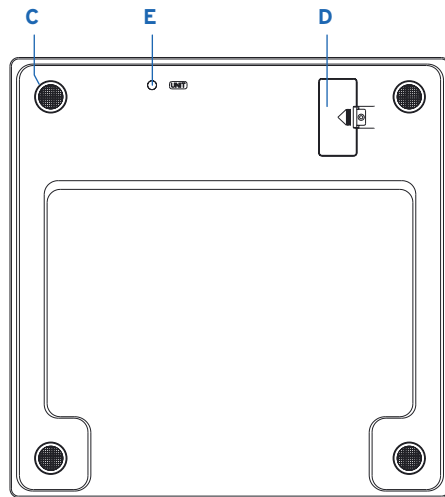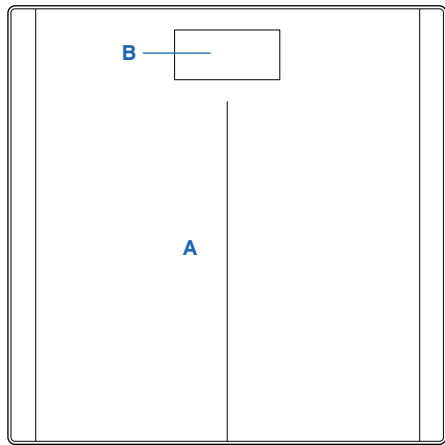Before First Use

1. Open the battery compartment on the underside of the scale. Remove the included AAA batteries. [\[Figure 1.1\]](#)
2. Remove the plastic strip.
3. Replace the batteries. Follow the guidelines in the battery compartment to match the + and - terminals. [\[Figure 1.2\]](#)

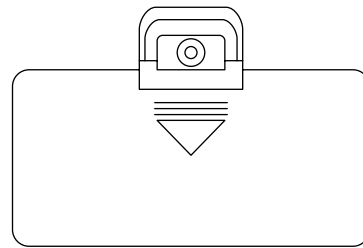
Figure 1.1

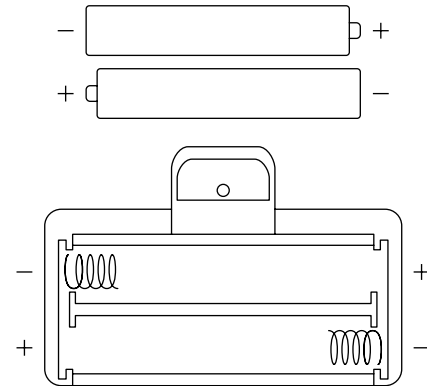
Figure 1.2

Operation

Using the Scale

Note: Press the unit switching button on the underside of the scale to switch between kilograms (kg), stone (st), and pounds (lb).

1. Place the scale on a hard, flat surface.

Note: Using the scale on soft surfaces such as carpets and mats will cause inaccurate measurements.

2. Step onto the scale. The scale will automatically turn on and start measuring your weight. [\[Figure 2.1\]](#)
3. Keep standing on the scale until the digits on the display flash 3 times, showing your final measurement. [\[Figure 2.2\]](#)

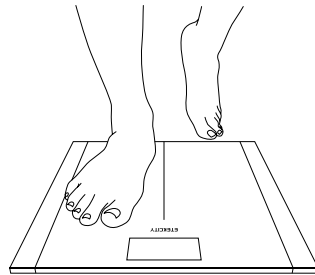
Figure 2.1

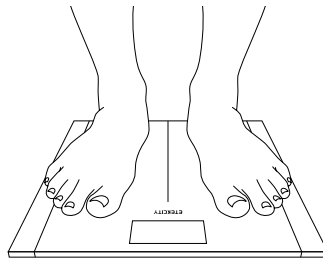
Figure 2.2

Body Measuring Tape

Use the included body measuring tape to track your fitness progress. Wrap the measuring tape around the body part you are measuring.

Replacing the Batteries

1. Open the battery compartment on the underside of the scale.
2. Remove the old batteries and dispose of them properly.
3. Install 2 new 1.5V AAA batteries. Make sure the polarities match. [\[Figure 3.1\]](#)
4. Close the battery compartment, place the scale on a hard, flat surface, and wait for the display to show "0.0" before using.


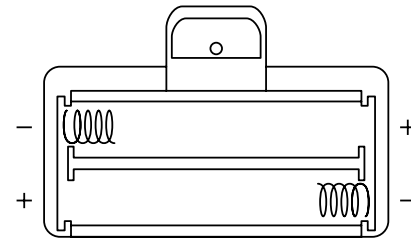
Figure 3.1

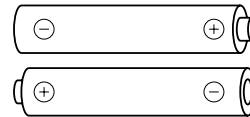Cleaning the Scale

1. Use a slightly damp cloth to clean the surface of the scale, followed by a soft, dry cloth to wipe away the remaining moisture.
2. **Do not** use chemical or abrasive cleaning agents to clean the scale.

Calibrating the Scale

If the scale has been moved or flipped upside down, it must be calibrated to ensure accurate results.

1. Place the scale on a hard, flat surface.
2. Step on the scale with one foot until the display lights up, then step off.
3. The scale will show "0.0", indicating that the calibration was successful.

Troubleshooting

Problem	Possible Solution
Scale does not turn on.	Open the battery compartment cover and remove the plastic strip.
	Replace the batteries (see page 13).
	Insert batteries under the correct polarity as marked in the battery compartment.
Display backlight does not come on when stepping on the scale.	Scale may be malfunctioning. Contact Customer Support (see page 17).
Display shows "Err".	Weight exceeds 180 kg / 400 lb. Use another scale with a higher weight range.
Display shows "LO".	Replace the batteries (see page 13).
Display shows "C".	Scale needs to be calibrated. Step off the scale, wait for the display to show "0.0", then step back on the scale.
Reading is inaccurate.	Place scale on a level surface.
	Step off the scale, wait for the display to show "0.0", then step back on the scale.
	Remove the batteries for 45 seconds and replace them to reset the scale.

If your problem is not listed, please contact [Customer Support](#) (see page 17).

Warranty Information

Terms & Policy

Etekcitec Corporation warrants all products to be of the highest quality in material, craftsmanship, and service for 2 years, effective from the date of purchase to the end of the warranty period. Warranty lengths may vary between product categories.

If you have any questions or concerns about your new product, please contact our helpful Customer Support Team.

Customer Support

Etekcitec Corporation

1202 N. Miller St., Suite A
Anaheim, CA 92806
USA

Email: support.eu@etekcity.com

**Please have your order invoice and order ID ready before contacting Customer Support.*

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