# Digital Body Weight Scale

**Model:** EB9380H Series

**User Manual**


EN

DE

FR

ES

IT
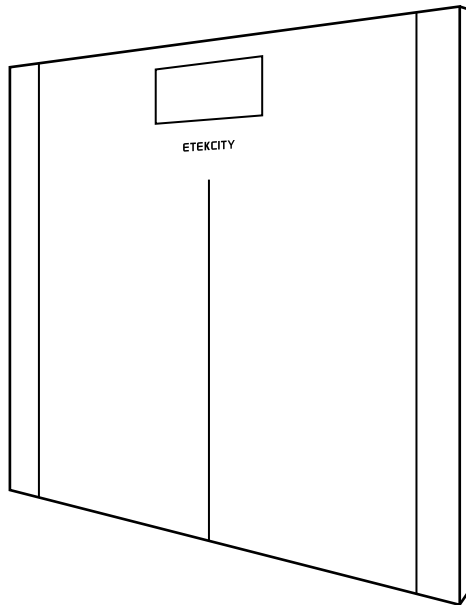
## Table of Contents

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**Thank you for purchasing the  
Digital Body Weight Scale by Etekcity.**

If you have any questions or concerns,  
please reach out to our helpful Customer  
Support Team at **support.eu@etekcity.  
com**. We hope you enjoy your new scale!

# Safety Information

**Please read and follow all instructions and safety guidelines in this manual.**

## General Use

- Check the scale before using. **Do not** use the scale if it is damaged in any way.
- **Do not** strike or bang on the scale and **do not** use the scale if there is a slight crack on the glass. Due to the nature of tempered glass, the glass will shatter if there is a slight crack. If you see a crack, contact **Customer Support** (see page 17).
- The glass platform is slippery when wet. **Always** make sure the platform and your feet are dry before stepping onto the scale.
- Place the scale on a hard, flat surface.
- To weigh safely and accurately, stand with your feet shoulder-width apart while on the scale. **Do not** stand on the edge of the scale or jump on the scale.
- Avoid standing on the scale for a long period of time.
- **Do not** overload the scale. Maximum weight capacity is 180 kg / 28 st / 400 lb.
- Keep the scale in a cool area with low humidity. Keep the scale away from heat sources, such as ovens or heaters.
- Store the scale horizontally, not vertically.
- For household use **only**. This scale is not intended for commercial use.
- Handle with care.



## Batteries

- **Do not** mix new and old batteries.
- If the scale is not in use for a long period of time, remove the batteries to extend the lifetime of the scale.

**SAVE THESE  
INSTRUCTIONS**

# Function Diagram

- A. Platform
- B. Display
- C. Anti-Skid Padding
- D. Battery Compartment
- E. Unit Switching Button
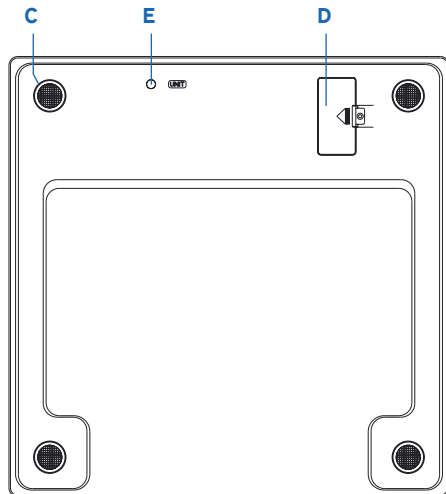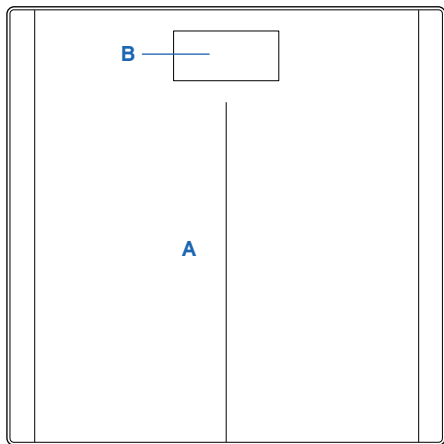
## Before First Use

1. Open the battery compartment on the underside of the scale. Remove the included AAA batteries. [Figure 1.1]
2. Remove the plastic strip.
3. Replace the batteries. Follow the guidelines in the battery compartment to match the + and - terminals. [Figure 1.2]

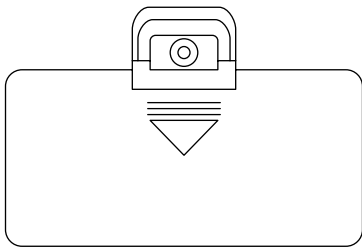
Figure 1.1

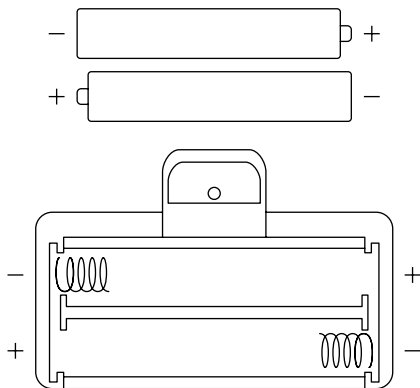
Figure 1.2

# Operation

## Using the Scale

**Note:** Press the unit switching button on the underside of the scale to switch between kilograms (kg), stone (st), and pounds (lb).

1. Place the scale on a hard, flat surface.

**Note:** Using the scale on soft surfaces such as carpets and mats will cause inaccurate measurements.

2. Step on the scale, then step off. Wait for the screen to display "0.0" before using.  
[Figure 2.1]
3. Stand on the scale until the display flashes several times, showing your final measurement.  
[Figure 2.2]

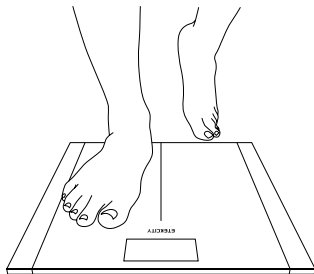
Figure 2.1

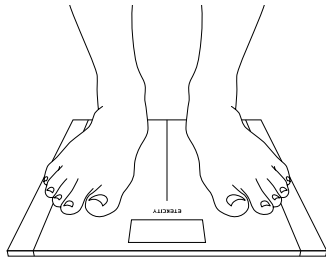
Figure 2.2

## Replacing the Batteries

1. Open the battery compartment on the underside of the scale.
2. Remove the old batteries and dispose of them properly.
3. Install 2 new 1.5V AAA batteries. Make sure the polarities match. [\[Figure 3.1\]](#)
4. Close the battery compartment, place the scale on a hard, flat surface, and wait for the display to show "0.0" before using.

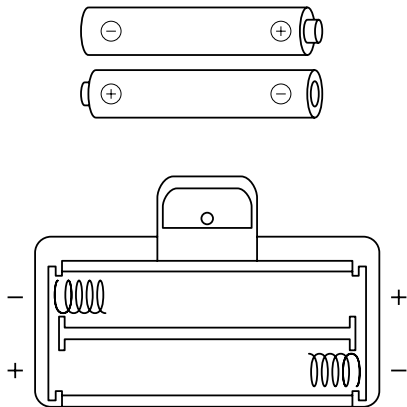
Figure 3.1

## Cleaning the Scale

1. Use a slightly damp cloth to clean the surface of the scale, followed by a soft, dry cloth to wipe away the remaining moisture.
2. **Do not** use chemical or abrasive cleaning agents to clean the scale.

## Clearing the Scale

If the scale has been moved or flipped upside down, it must be cleared to ensure accurate results.

1. Place the scale on a hard, flat surface.
2. Step on the scale until digits appear on the display, then step off.
3. The scale will show "0.0" when cleared.

# Troubleshooting

Problem	Possible Solution
<b>Scale does not turn on.</b>	Open the battery compartment cover and remove the plastic strip.
	Replace the batteries (see page 13).
	Insert batteries under the correct polarity as marked in the battery compartment.
<b>Display backlight does not come on when stepping on the scale.</b>	Scale may be malfunctioning. Contact <a href="#">Customer Support</a> (see page 17).
<b>Display shows "Err".</b>	Weight exceeds 180 kg / 400 lb. Use another scale with a higher weight range.
<b>Display shows "LO".</b>	Replace the batteries (see page 13).
<b>Display shows "E".</b>	Scale needs to be cleared. Please place the scale on a flat surface and wait for it to be cleared before using it.
<b>Reading is inaccurate.</b>	Place scale on a level surface.
	Step off the scale, wait for the display to show "0.0", then step back on the scale.
	Remove the batteries for 45 seconds and replace them to reset the scale.

If your problem is not listed, please contact [Customer Support](#) (see page 17).

# Warranty Information

## Terms & Policy

Etecity Corporation warrants all products to be of the highest quality in material, craftsmanship, and service for 2 years, effective from the date of purchase to the end of the warranty period. Warranty lengths may vary between product categories.

If you have any questions or concerns about your new product, please contact our helpful Customer Support Team.

# Customer Support

## **Etekcitec Corporation**

1202 N. Miller St., Suite A  
Anaheim, CA 92806  
USA

**Email:** [support.eu@etekcity.com](mailto:support.eu@etekcity.com)

*\*Please have your order invoice and order ID ready before contacting Customer Support.*

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## **Etekciti Corporation**