

QUICK START GUIDE

SPECIFICATIONS

Weight Capacity	400 lb / 180 kg
Increments	0.2 lb / 0.1 kg
Weight Units	lb / kg
Platform	Tempered glass
Dimensions	11.8 x 11.8 x 0.95 in / 30 x 30 x 2.4 cm
Battery	2 x 1.5V AAA batteries
Auto-Off	15 seconds

PACKAGE CONTENTS

- 1 x Digital Body Weight Scale
- 2 x 1.5V AAA Batteries (Pre-Installed)
- 1 x Measuring Tape
- 1 x Quick Start Guide

TIPS


- To switch between pounds (lb) and kilograms (kg), press the unit switching button on the underside of the scale.
- If your scale has been moved or flipped upside down, it will need to be calibrated.

USING YOUR SCALE

Note: Before first use, remove the plastic strip from the battery compartment.

1. Place scale on a hard, flat surface (not a carpet or mat).
2. Step onto the scale. The scale will automatically turn on and start measuring your weight.
3. Stand on the scale until the digits on the display flash 3 times, showing your final measurement.

1. Step on the scale until digits appear on the display, then step off.
2. The scale will show “**[**” while calibrating, then “**0.0**” when complete.

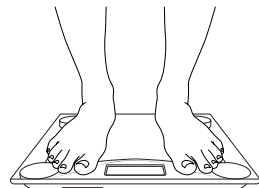

ONLINE MANUAL



Need more help? Scan here to view the manual for more information.

You can also type the following link in a web browser:
www.etekcity.com/support/manual/model/eb4074c

WARRANTY



Etekcity Corporation warrants all products to be of the highest quality in material, craftsmanship, and service, effective from the date of purchase to the end of the warranty period. Warranty lengths may vary between product categories.

Scan the QR code to register your product and extend your 1-year warranty by an additional year.

CUSTOMER SUPPORT

Etekcity Corporation

1202 N. Miller St., Suite A
Anaheim, CA 92806

Email: support@etekcity.com

Toll-Free: (855) 686-3835

Support Hours

Monday-Friday
9:00 am-5:00 pm PST/PDT

*Please have your invoice and order ID ready before contacting Customer Support.